

# MARSDEN

M-650

M-651

M-652

M-653



Class III/MDD Approved platforms for weighing wheelchairs

Fast, simple wheelchair weighing accurate to 100g

Available with handrails or indicator column

[www.marsden-weighing.co.uk](http://www.marsden-weighing.co.uk)



# Weighing a patient in a wheelchair shouldn't be more difficult than weighing one that can stand.

It's often taken as a given that weighing a patient unable to stand is more complicated and time consuming than weighing one that can stand.

But why is this? As long as the scale itself is accessible enough, there should be no reason to struggle weighing a patient in a wheelchair. **The M-65X series has been designed so that there is no struggle.**

## Marsden M-65X series: Four solutions in one scale

Because every patient is different, and every weighing requirement is different, we have created **four solutions from one weighing scale.**

The M-650 is the simplest of our wheelchair scales. It is a large flat area, with **shallow ramps on two sides**, that allows almost any wheelchair to be rolled onto it and weighed. **Class III Approved 100g accuracy** ensures it's capable of providing **weight readings for almost any medical need.**

Its **easy to use indicator** makes getting an accurate weight reading quick and simple. The **Tare** function means you can remove the weight of the wheelchair from the reading. Enter the patient's height using the **full numeric keypad** and you'll get their **BMI** reading, too.

Then, when you've finished using the M-650, simply lift it onto its side using the handle, and roll away using the integrated wheels.

## Indicator at waist height for easier operation

For extra convenience, the M-651 has all the features of the M-650 but has a **column-mounted indicator.**

This makes operation of the scale easier, particularly if you're intending to use it in a more fixed position.



## Do patients want to stand? Choose handrails...

Fitted handrails make our wheelchair scale perfect for those who wish to stand, but need a bit of support.

Choose **two handrails** (M-652) or go for the M-653, with a **single folding handrail**, perfect for transportation and storage. The indicator is conveniently handrail-mounted.

## Key features

300kg capacity; accuracy to 100g

MDD Approved

Class III Approved

Powered by rechargeable battery

Hold/Tare/BMI

Low profile base with shallow ramps

Handle and wheels for portability

Column-mounted indicator (M-651)

Two handrails (M-652)

One hinged handrail (M-653)

Platform dimensions: 1150 x 800 x 66mm

## Record weights or send them to patient records

So, the M-65X series is a versatile range of solutions for your weighing needs. But what they offer doesn't stop there.

Each of the four variants feature our **premium indicator**, which is available with **Bluetooth, Wifi or Micro USB connectivity**. Choose one of these options and, with your patient database configured correctly, **weight, height and BMI data can be populated instantly in their records**.

Plus, you can choose the **optional printer** and get an **instant print-out of the patient's readings**.



WEIGHT	97.45kg
TARE WEIGHT	0.00kg
NET WEIGHT	97.45kg
PATIENT HEIGHT	178.8cm
PATIENT B.M.I	30.4
14/08/2014	10:52

 **MARSDEN**  
ACCURACY ASSURED

Print

# Key features at-a-glance

Column-mounted indicator

Wheels and handle  
(common to all models)

Low profile base with shallow ramps  
(common to all models)

Multi-functional indicator with Tare, Hold, BMI  
(common to all models)

One hinged handrail


Two handrails



M-652

M-653

# Weighing area dimensions

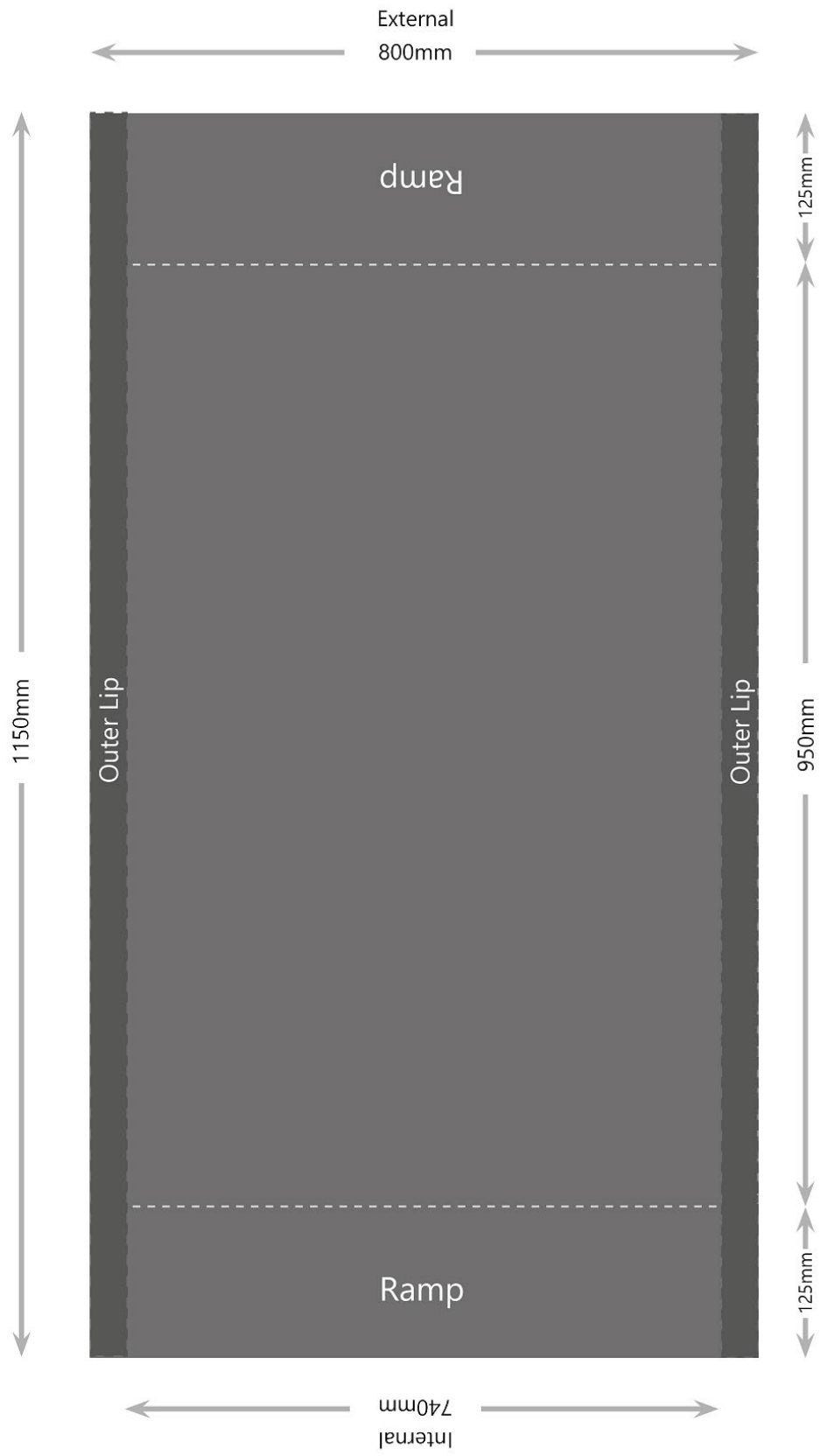

# MARSDEN

Tel: 01709 364296 / 0800 169 2775

Fax: 01709 364293

E-mail: [sales@marsdengroup.co.uk](mailto:sales@marsdengroup.co.uk)

[www.marsden-weighing.co.uk](http://www.marsden-weighing.co.uk)

 MarsdenWeighing

 @MarsdenWeighing

 [www.youtube.com/MarsdenWeighing](http://www.youtube.com/MarsdenWeighing)

 <http://pinterest.com/marsdenweighing>