


MARSDEN

H-629-2 Body Measuring Rod

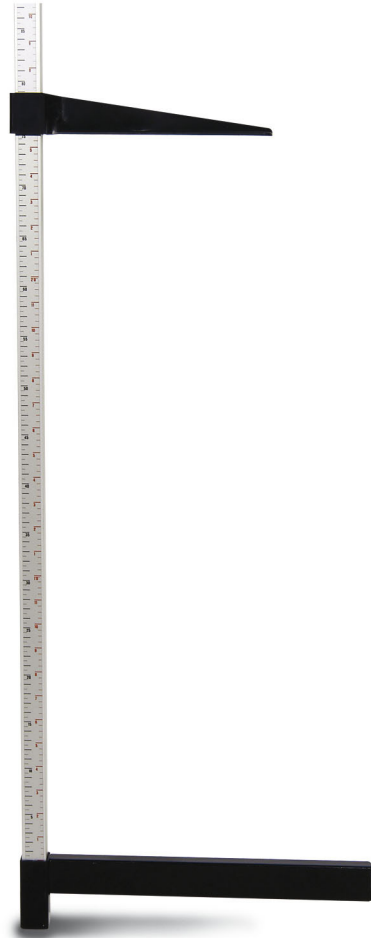

- Lightweight aluminium construction
- Sliding cursor for accurate measuring
- 1mm increments
- Plastic foot support

Marsden H-629-2: Body Measuring Rod

- Lightweight aluminium construction
- 100cm measuring range
- Sliding cursor for accurate measuring
- 1mm increments
- Fitted with a foot support
- Ideal for use in mortuaries
- 1m version also available
- 1 year warranty as standard



Multi-use height measure

The H-629-2 is a multi-functional body measuring rod. This device is designed to be held against a person whilst they are in a horizontal position.

Since the release of the H-629 product range we have found them to be increasingly popular in hospital wards with many immobile patients, baby clinics and within mortuaries.

Foot support and cursor

Due to our diligent research and development upon release of this height measure, we added a foot support as standard.

We found that when measuring people in a horizontal position, the height measure had to have an anchor point to ensure that an accurate height was obtained and the measure did not move.

Easy to use

The H-629-2 is incredibly simple and easy to use, with only one moving part. Simply ensure that the patient is laying perfectly flat and their bed is not tilted.

Place the height measure parallel to the patient, with the foot support resting on the base of the patients feet, then adjust the cursor until it is sitting flat on top of the head.

Measuring Range: 0cm to 200cm

Dimensions: 2000mm to 350mm

Graduations: 1mm

Weight: 0.6kg