


# MARSDEN

## M-545 Entry Level Floor Scale with BMI

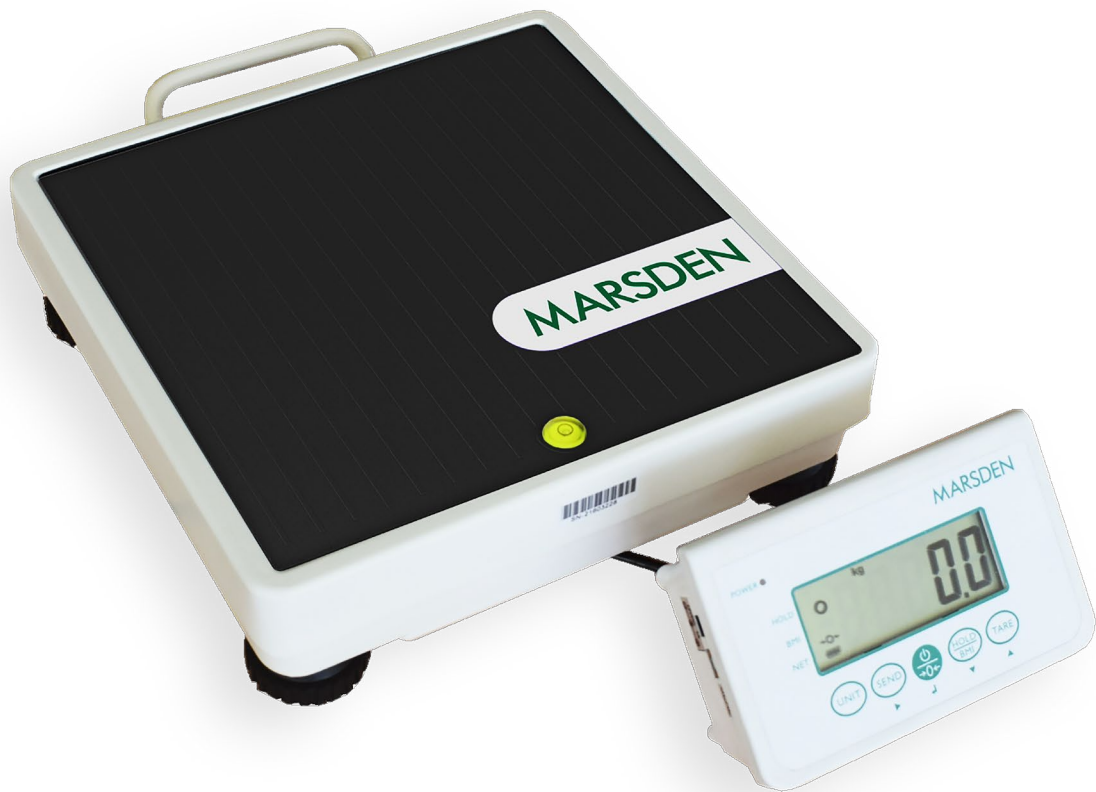

- Class III and MDD Approved
- Accurate to 100g
- Tough, robust and hard wearing base
- Built in carry handle for portability
- Suitable for GPs and hospitals
- Optional carry case available

# Marsden M-545:

## Entry Level Floor Scale

- Class III and MDD Approved
- Accurate to 100g
- Perfect for use in the GP surgery or hospital
- Powered by 6 x AA batteries
- 160 hours of continuous use
- Auto power-off function helps to save battery power
- Simple and easy to use
- Optional carry case available to aid portability
- Easy to read remote display
- Hold and Tare functions to enable accurate and easy weighing
- Made from reinforced aluminium structure
- Covered by the Marsden 8 year warranty



### Simple weighing

The M-545 is a simple and easy to use stand on scale, made from reinforced aluminium, which is ideal for weighing all patients up to a capacity of 250kg.

The scale has been designed to be hard wearing making it ideal for heavy use in busy traffic areas, such as general wards, outpatient departments and busy GP surgeries.

### Remote display as standard

The M-545 comes with a remote display fitted as standard. This means that when you are operating the scale in sensitive environments, the weight reading is discreet.

The remote display means that the user of the scale can see the weight reading with ease and enabled the scale to be fully operable without the need to bend down and press buttons.

### Easy to use

The M-545 features our vastly popular I-300 indicator. This has been a best-seller for us due to the 5 simple buttons and the bright, easy to read display.

The On/Off button ensures that you have full control over the device. Hold enables you to gain weights even for the patient is not steady on the scale and the tare key allows you to deduct any unwanted weight.

Capacity:	250kg	Power Supply:	Mains/6 x AA Batteries
Graduations:	100g	Battery Life:	160 hours of continuous use
Functions:	Hold, Portable, Tare	Approvals:	Class III and MDD Approved
Dimensions:	360mm x 310mm x 90mm	Optional Extras:	Carry Case
Weight of Scale:	6kg	Warranty Period:	8 Years