

# MARSDEN

## M-565 High Capacity Slimming Scale




- Suitable for use in gyms and slimming clubs
- Accurate to 200g/0.5lb
- Tough, robust and hard wearing base
- Built in carry handle for portability
- USB Type B port to enable data transfer
- Optional carry case available

# Marsden M-565:

## Entry Level Floor Scale

- Non-approved to enable weighing in stones and pounds
- Accurate to 200g/0.5lb
- Perfect for use in gyms and health clubs
- Powered by 6 x AA batteries
- 180 hours of continuous use
- Auto power-off function helps to save battery power
- Simple and easy to use
- Optional carry case available to aid portability
- Easy to read remote display
- Hold and Tare functions to enable accurate and easy weighing
- Made from reinforced aluminium structure
- Covered by the Marsden 4 year warranty

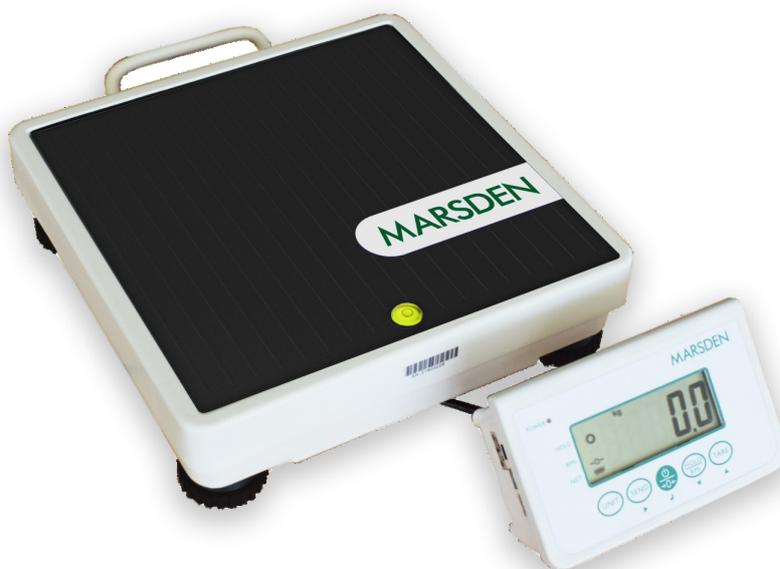

### Simple weighing

The M-565 is a simple and easy to use stand on scale, made from reinforced aluminium, which is ideal for weighing all patients up to a capacity of 250kg/39st.

The scale has been designed to be hard wearing making it ideal for heavy use in busy traffic areas, such as general wards, outpatient departments and busy GP surgeries.

### Remote display as standard

The M-565 comes with a remote display. This means that when you are operating the scale in sensitive environments, the weight reading is discreet for only the user to see.

The remote display also means that the user does not have to bend down and see the weight reading on the display, as it is easy to see on the device. Plus, the device is entirely controlled via the remote display.

### Easy to use

The M-565 features our popular I-300 indicator, which is very simple to use due to the 5 buttons and a large clear display which is easily visible.

The On/Off button ensures that you have full control over the device. Hold key enables you to gain weights even for the patient is not steady on the scale and the tare key allows you to deduct any unwanted weight.

Capacity:	250kg/39st	Power Supply:	6 x AA Batteries
Graduations:	100g/0.5lb	Battery Life:	180 hours of continuous use
Functions:	Hold, Portable, Tare	Approvals:	Class III and MDD Approved
Dimensions:	360mm x 310mm x 90mm	Optional Extras:	Carry Case
Weight of Scale:	6kg	Warranty Period:	1 Year