

MARSDEN

M-615 Lightweight Folding Wheelchair Scale



- Class III Approved
- Accurate to 100g
- Portable, folding wheelchair beams
- Ideal for use in the community
- Optional carry case with wheels available
- Optional Wi-Fi or Bluetooth available

Developed in partnership with:




**Birmingham
Community Healthcare**
NHS Foundation Trust

Marsden M-615:

Portable, Folding Wheelchair Weigh Beams

- Class III Approved
- Accurate to 100g
- Perfect for use in the community and hospital
- Powered by rechargeable battery pack
- 55 hours of continuous use
- Auto power-off function helps to save battery power
- Simple and easy to use
- Easy to read display
- Hold and Tare functions to enable accurate and easy weighing
- Optional carry case with wheels available
- Designed with input from Birmingham Community Healthcare NHS Foundation Trust

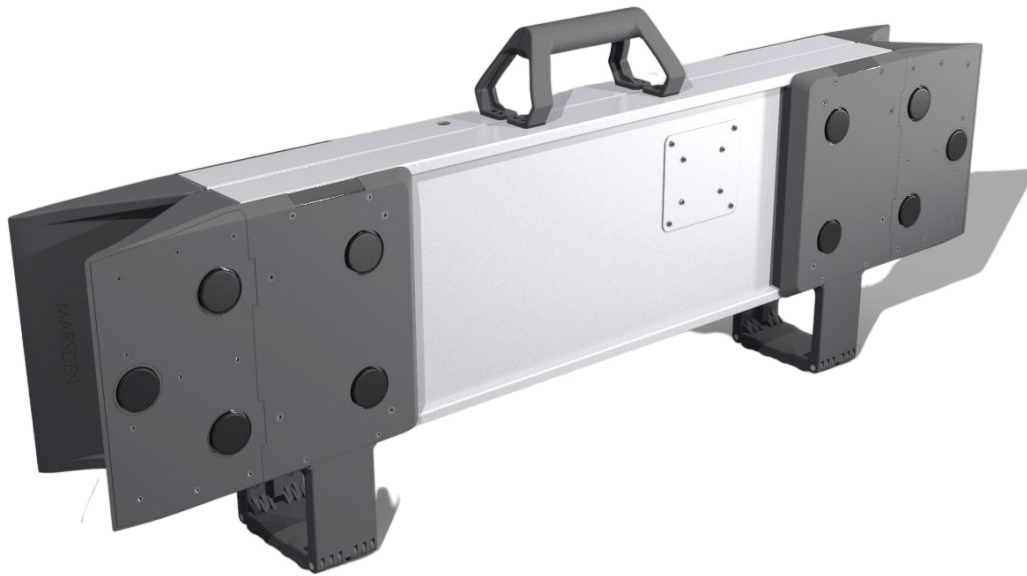

Folding design

The M-615 was designed in collaboration with Birmingham Community Healthcare NHS Foundation Trust to ensure that accurate weights can be easily obtained.

The folding design of the M-615 weigh beams means that the scale is easy to transport and incredibly easy for the people to set up. Simply unfold the beams, connect the indicator and you are ready to weigh.

Optional connectivity options

The M-615 is available with either Bluetooth or Wi-Fi connectivity to ensure that our customers are able to transfer data seamlessly or input into EPR systems.

We would simply need to insert a connectivity module at the time of order and then provide you with our standard protocols for easy and simple integration into your existing system.

Easy to use

The M-615 features from our best-selling DP-3810 indicator. This device features from a full numeric keypad for easy BMI calculation and simple function keys.

The On/Off button ensures that you have full control over the device and you can prolong the battery life. Hold enables you to gain readings for unsteady patients and Tare allows you to deduct unwanted weight.

Capacity:	300kg	Power Supply:	Mains/Rechargeable
Graduations:	100g	Battery Life:	55 hours of continuous use
Functions:	Hold, Portable, Tare	Approvals:	Class III Approved
Dimensions:	200mm x 740mm	Optional Extras:	Bluetooth/Wi-Fi, Carry Case, Printer
Weight of Scale:	9.8kg	Warranty Period:	4 Years