


MARSDEN

M-930 Pit Mounted Bed and Trolley Weigher

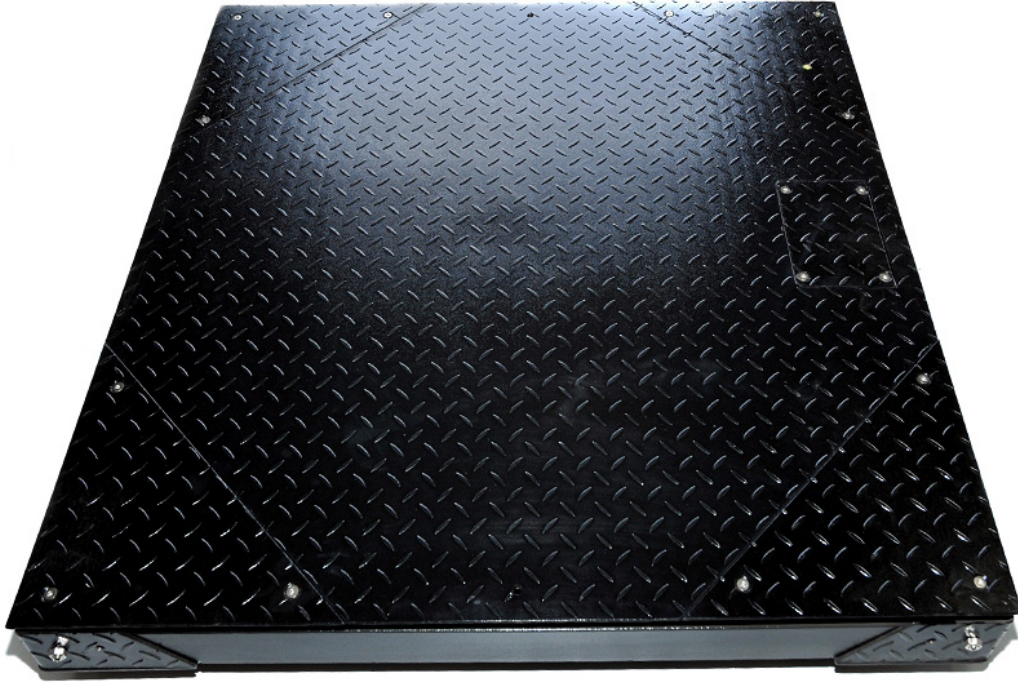


- Class III Approved
- Accurate from 200g
- Heavy duty weighing platform
- Suitable for weighing beds in hospitals
- Level with the floor for seamless weighing
- Optional Wi-Fi or Bluetooth available

Marsden M-930:

Pit Mounted Floor Scale

- Class III Approved
- Accurate from 200g
- Perfect for weighing beds and trolleys in hospitals
- Powered by rechargeable battery pack
- 55 hours of continuous use
- Auto power-off function helps to save battery power
- Flush with the floor for seamless weighing
- The patient will not feel that they are on a weighing scale
- Easy to read display
- Hold and Tare functions to enable accurate and easy weighing
- Optional Bluetooth and Wi-Fi available
- Covered by the Marsden 4 year warranty

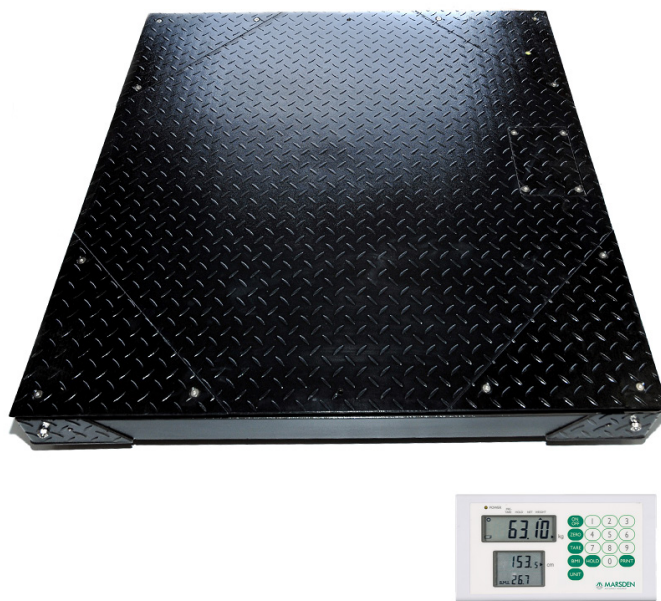

Why the M-930?

The M-930 is a pit mounted bed weighing scale, perfect for weighing bedridden patients or those in trolleys upon entrance to the hospital.

We will work closely with your internal hospital departments to ensure that any preliminary work is carried out effectively and to our specifications. Marsden engineers will then carry out installation.

Optional connectivity options

The M-930 is available with either Bluetooth or Wi-Fi connectivity to ensure that our customers are able to transfer data seamlessly or input into EPR systems.

We would simply need to insert a connectivity module at the time of order and then provide you with our standard protocols for easy and simple integration into your existing system.

Easy to use

The M-930 features from our best-selling DP-3810 indicator. This device features from a full numeric keypad for easy BMI calculation and simple function keys.

The On/Off button ensures that you have full control over the device and you can prolong the battery life. Hold enables you to gain readings for unsteady patients and Tare allows you to deduct unwanted weight.

Capacity:	600kg or 1500kg	Power Supply:	Mains/Rechargeable
Graduations:	200g or 500g	Battery Life:	55 hours of continuous use
Functions:	Hold, Tare, BMI	Approvals:	Class III Approved
Dimensions:	1200mm x 1200mm / 1500mm x 1500mm	Warranty Period:	4 Years